


Certificate of Analysis
No. CA22169

Product: Eichrom's UTEVA® Resin UT15-C20-S-3PA Columns
 Lot #: UT22E
 Date of Manufacture: 21-JUL-2022

Eichrom's UTEVA® Resin UT15-C20-S-3PA column Lot UT22E has been manufactured according to Eichrom's standard operating procedures. The resin used to manufacture this lot, UTS041222, has been tested in accordance with Eichrom's quality assurance policy and procedures. This lot of product meets all Eichrom quality specifications.

ANALYSIS RESULTS

<u>Test Value</u>	<u>Specification</u>	<u>Result</u>
Appearance	White beads	Pass
Distribution Ratio, 3.1 N nitric acid	>100	165.7
% ²³⁰ Th stripped in 5M HCl/0.05M oxalic acid		
Least acceptable mean (4 rep)	90.0	99.1
Max. standard deviation	10.0	1.5
% ²³⁰ Th stripped in 1M HCl		
Maximum acceptable mean (5 rep)	0.5	0.1
Max. standard deviation	1.0	0.0
% ²³³ U stripped in 1M HCl		
Least acceptable mean (5 rep)	90.0	97.7
Max. standard deviation	10.0	1.6
% ²³³ U stripped in 5M HCl/0.05M oxalic acid		
Maximum acceptable mean (5 rep)	1.0	0.4
Max. standard deviation	1.0	0.2

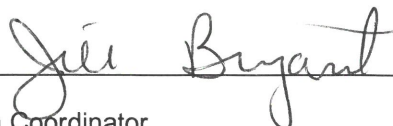

 Jill Bryant
 Quality System Coordinator

July 21, 2022